

FREEDOM 5000

SPECIFICATIONS



SEAT WIDTH
38cm, 40.5cm, 43cm,
45.5cm, 48cm, 50.5cm



SEAT DEPTH
42.5cm



OVERALL WIDTH
Seat Width - 20 cms



SEAT TO FLOOR HEIGHT
46cm - 51cm



BACKREST ANGLE
7 degrees



BACKREST HEIGHT
43 cms



OVERALL LENGTH
(incl. footrest)
103cm - 109cm
(Self Propel)
95cm (Transit).



WEIGHT CARRYING CAPACITY
125 kgs



NET WEIGHT
17.7 kg (approx)



GROSS WEIGHT
19.7 kg

FRAME COLOURS



Anthracite Black Pearl Grey Ocean Blue

- All plastic parts are made with Nylon PA66-30% GF & PA6-30% GF for better durability.
- It comes with a 2 year warranty on the main frame, Crossbars and other detachable parts.



Registration, Evaluation, Authorisation and Restriction of Chemicals



Restriction of Hazardous Substances



Certified Flame Retardant Fabric



UK approved crash testing

FREEDOM
RESTORE LIFE



FREEDOM 5000

VARIANT: SELF PROPELLED

LIGHTWEIGHT | STEEL BODY | FOLDABLE


FREEDOM 5000

is the ideal solution for those looking for a quality, robust, economical and easy-to-use manual wheelchair.

FREEDOM 5000 is a strong and lightweight foldable steel wheelchair that has been designed for easy transportation. The frame is easily adjustable that makes it suitable for daily use. A safety belt and heel strap have been incorporated into the design for safety and the seat comes with an attachable padded cushion for extra comfort, making this model a perfect aide to face any challenge head-on.

VARIANTS AVAILABLE

1 SELF PROPELLED WITH 24" REAR WHEEL

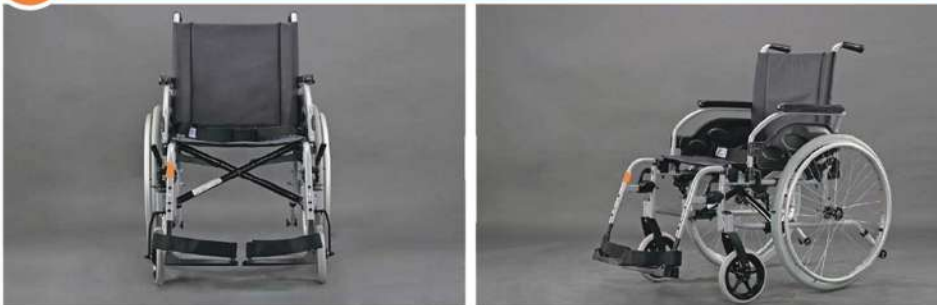

2 TRANSIT WITH 12" REAR WHEEL



FREEDOM 5000

FEATURES



DETACHABLE, SWING IN SWING OUT LEGREST



HEIGHT ADJUSTABLE FOOT TUBE



EFFORTLESS AND CONVENIENT FOLDING



BRAKE LEVER HEIGHT ADJUSTMENT FOR EASY USER SHIFTING



DETACHABLE, SWING AWAY ARMRESTS



DETACHABLE 24" PU TYRE WITH QUICK RELEASE AXLE



ANTI TIPPERS FOR SLOPES



PADDED CUSHION



VARIANT: TRANSIT

YOUR DAILY COMPANION

- Steel, foldable, with aluminum forks.
- Tricolour chassis (black, blue and gray).
- Cross braces (X-members) for strut support between the two side frames.
- Compact and foldable for easy transport and storage.
- Padded cushions and nylon backrest for ultimate comfort.
- Backrest fixed at 7° for optimum ease.
- Foldable armrests with upholstered cushion.
- Detachable foldable footrests for ease-of-use.
- 12" or 24" rear wheels (hand rim in anodized aluminum).
- 8" front wheels.
- Available with a wide range of accessories.