

SSS100

Self Propelled Shower
and Commode Wheelchair



Product Summary (SSS100)

The Arcatron SSS100 is a self propelled wheelchair designed for senior citizens and individuals with limited mobility to access bathrooms in a dignified and safe manner. Our unique tool free height adjustability eases transfers to bed and allows rolling in over different commode heights. Easy to disinfect, waterproof cushions enable to maintain hygiene standards. The big wheels can be detached with a click of a button with quick release axle. SSS100 frame is corrosion resistant and designed for a comfortable and safe experience.

1. Made of Non Rusting Materials



2. Tool free height adjustability



3. Front seat cutout for easy hygiene access



4B. Chair with top seat cover*



4A. Top Seat Cover

Can be used as a regular wheelchair with the top seat cover

* top seat cover is an optional accessory



5. Tool Free Assembly

Five Levels of Tool Free Height Adjustment



Rolls Over Commodes

Easy Bed Transfer with Height Adjustment and Swivel Armrests



'One Push' - Add or Remove Big Wheels Instantly!
Mag wheels with quick release axle.

Can be used as an
Attended Propelled Chair.
(SAS100)




Mag Wheels can be
inserted instantly at
5 different heights

SSS100

SPECIFICATIONS (Self-Propelled)

PARAMETERS	NOTATION	MODEL	
		SSS100	SSS100N (Without Magwheel Pushrims)
Dimension	H* x W x D	34-38" x 26" x 42"	34-38" x 24-25.5" x 42"
Weight	-	22.75 Kg	22.75 Kg
Weight Capacity	-	120 Kg	120 Kg
Seat Dimensions	C x A	19" x 16"	19" x 16"
Width Between Arms	E	16.5"	16.5"
Commode Clearance	M* x G	(17"-21") x 16.5"	(17"-21") x 16.5"
Seat Depth	A	16"	14.5-16"
Foot Plate Distance	J	19-21"	19-21"
Seat to Arm Height	F	8"	8"
Seat Height	B*	20" - 24"	20" - 24"
Back Dimensions	K x L	17" x 16"	17" x 16"

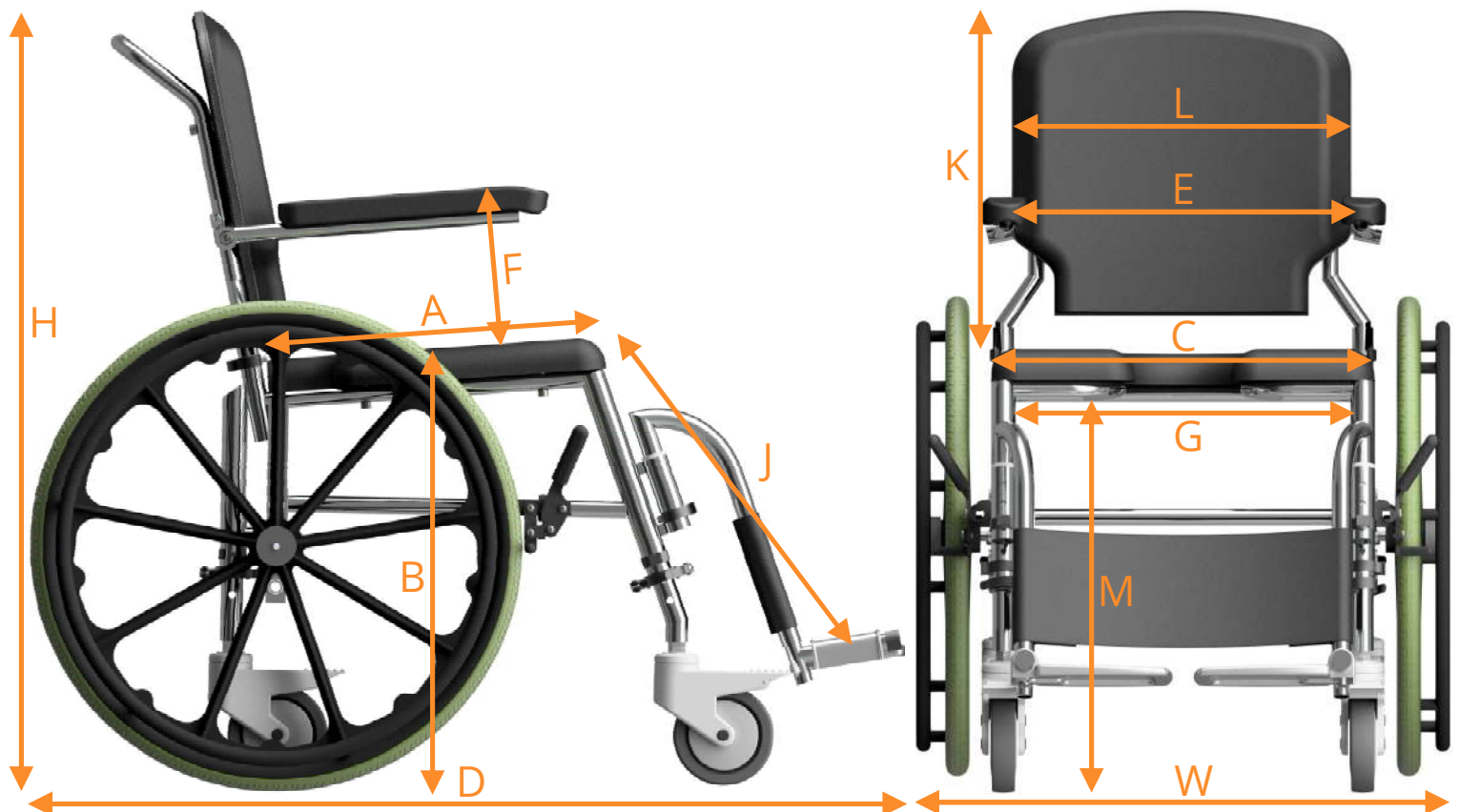

*For 19" Commode clearance



+91 9049291555
+91 20 66106643

contact@arcatron.com

www.arcatron.com